

# OPEN CARE STATION


SPN-1600 / SPN-1900 / SPN-1900Z



UNIMED

## OPEN CARE STATIONS SPN-1600 / SPN-1900 / SPN-1900Z

Infant Care Stations type SPN-1600, SPN-1900 and SPN-1900Z are designed to be used for nursing and care treatments of infants in „rooming in” rooms and traditional neonatal hospital units. It provides integration of all actions related to nursing newborn in one place. Every station is equipped with nursing table, bathtub, sink and infant warmer (infrared radiator emitting infrared radiation) type IR-10, which maintain proper temperature parameters during care treatments. Also SPN-1900Z is equipped with sink and rotating faucet with shower. The SPN stations allows to perform efficiently and effectively following nursing procedures: warming, examining, weighting bathing, changing diapers and clothes, other newborn cosmetology, and other routine treatments. In option, infant care station can be equipped with newborn scales.

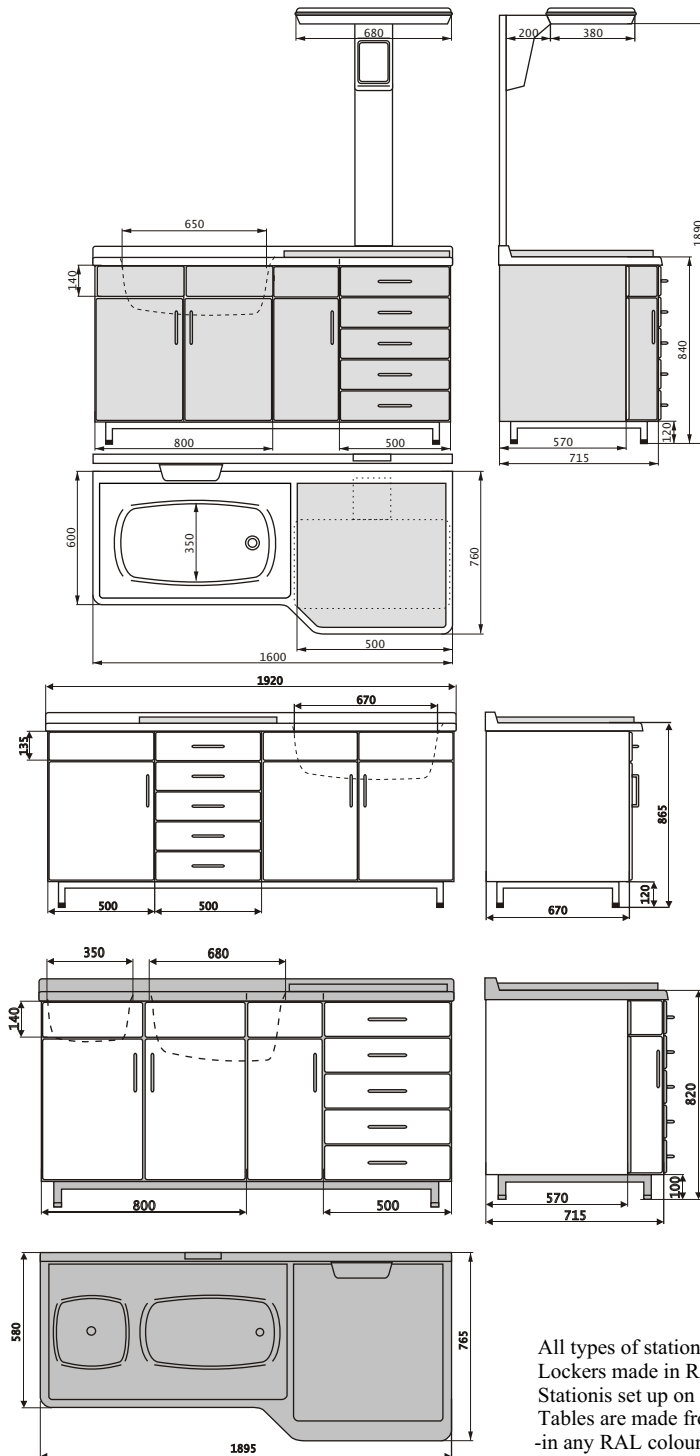

### OPEN CARE STATION SPN-1600

Dimensions of station shown on sketch

- Refuse bin with automatically lifting cover
- Lockers with shelves
- PCV covered mattress in colour matched with table
- Table with bathtub made of glassfiber -polyester laminate in any RAL colour scale
- Faucets (standard, optional with shower and thermostat)
- Solid base on heels
- Infant warmer type IR-10 emitting Infrared (optional fully automatic air temperature regulation) with body temperature sensor
- Column fixing warmer to station
- Optional drawer for scales.

Open care Station SPN-1600 are made in two versions SPN-1600P with bathtub on the right side of the station and SPN-1600L with bathtub on the left side of the station

### OPEN CARE STATION SPN-1900

Dimensions of station shown on sketch

- Locker with 3 pull-out basket drawers
- Pull-out plastic drawers 5 pcs.
- Refuse bin with automatically lifting cover
- Locker with shelves
- PCV covered mattress in colour matched with table
- Table with bathtub made of glassfiber -polyester laminate in any RAL colour scale
- Faucets (standard, optional with shower and thermostat)
- Solid base on stands
- Infant warmer type IR-10 emitting Infrared (optional fully automatic air temperature regulation) with body temperature sensor
- Column fixing warmer to station
- Place to put electric scales.

### OPEN CARE STATION SPN-1900Z

Dimensions of station shown on sketch

SPN-1900Z is a version of SPN-1900 equipped with additional sink and rotating faucet with shower. It is designed for rooms without stationary washbasins.

All types of stations are equipped with MDF lockers and drawers with round edges. Lockers made in RAL 9016 colour. Stations set up on metal frame with regulated heels. Tables are made from polyester laminate strengthened by glass fiber, covered with sanitary gelcoat -in any RAL colour scale. All stations are equipped with infra red radiant infant warmers IR-10 type and optional electrical scales Radwag WPT 6/15D.